



CX ORTHO INSTANT TREATMENT PEARLS!

In future mailings we'll periodically send a quick pearl or two that I've utilized in my practice. Hopefully these Pearls will save you and your practice some time and money!

.012 NITI and Bonded Eyelets

I've discussed before how much I LOVE .012 NITI! I use it 99.9% of the time as my initial arch wire. It aligns teeth super-fast and I feel like it is less painful for the patient as well.

Combined with bondable eyelets, I've seen some amazing results soon after the initial bonding!

Case AM Pre Treatment



Case AM After 16 weeks (2 visits)



Note on the above case that Large Opening Eyelets were used on Lower 4's, Lower Left 2, and Small Opening Eyelet was used on Lower Right 1. Also, the 2nd Photo shows a .016 NITI that was just placed prior to taking photo.

Things to be aware of with .012 as your initial arch wire:

- I rarely engage the 6's as the wire is so flexible it almost guarantees an emergency wire out.
- If patients chew ice the wire can easily freeze and fracture.

VISIT STORE

Need eyelets?

[Click Here](#)

Need .012 Niti wire?

[Click Here](#)




\$6.95 Priority Shipping

We ship same day when orders are placed before Noon Pacific Time.



Satisfied or refunded

We offer full refunds for any reason if not satisfied.

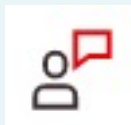

Support

Call us or email us and we'll answer your questions promptly.



Secure payments

Your payment info is safe and protected.



[Click here](#) to unsubscribes or manage your preferences

© Cx Orthodontic Supply. All Rights Reserved

961 E Nelson Rd Moses Lake, WA 98837, United States.

GO TO STORE >