



July Archwire Sale

10% Off on All Archwire, Offer Ends On 31-July

Have you tried using .012 NITINOL as your initial archwire?

We've used .012 Nitinol as our first wire in our practice for over 25 years with great success. Teeth align super-fast, and patients seem to report less pain with this initial wire.

Here are some tips when starting with .012 Nitinol:

1. At the full bonding appointment, we only tie in the archwires U/L 5-5. This avoids emergency appointments from the wire pulling out of the 6's due to the long inter bracket distance in those areas.
2. Reinforce to patients not to chew gum or ice with this flexible wire to avoid breakage.
3. On the 2nd appointment, we often ligature tie any remaining rotations. The flexible wire allows full engagement of the rotated tooth.
4. When using .012 for partial bonding (ie-U 2-2) we always crimp or bend the end of the wire to avoid poking the patient's lip or cheeks as the teeth align.



Because we love .012 Nitinol so much we are having a 10% sale on ALL ARCHWIRES until July 31!

VISIT STORE

[Click Here](#)

to check out our ARCHWIRES and get 10% off!




\$6.95 Priority Shipping

We ship same day when orders are placed before Noon Pacific Time.



Satisfied or refunded

We offer full refunds for any reason if not satisfied.

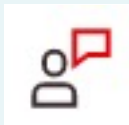

Support

Call us or email us and we'll answer your questions promptly.



Secure payments

Your payment info is safe and protected.



[Click here](#) to unsubscribes or manage your preferences

© Cx Orthodontic Supply. All Rights Reserved

961 E Nelson Rd Moses Lake, WA 98837, United States.

GO TO STORE >